

ENGINEERING PHARMACEUTICAL INNOVATION



Centrifugal and Positive Displacement Pumps for Hygienic Applications

ISPE Meeting November 18, 2008

Presented by:

Allen E. LeBoeuf, Product Manager

Oliver M. Dean, Inc. – Worcester, MA

Primary Issues:

- Selection of Technology.
- Control of Pumps.
- Materials of Construction / Certifications.
- Mechanical Seals.
- Suction Conditions.
- Efficiency / Power Consumption.
- Cleanability / Drainability / CIP-ability.

ENGINEERING PHARMACEUTICAL INNOVATION



Selection of Technology

When to Apply a Centrifugal Pump.

When to Apply a Positive Displacement .
(Discussion limited to rotary lobe style)

ENGINEERING PHARMACEUTICAL INNOVATION



Centrifugal vs. Positive Displacement



ENGINEERING PHARMACEUTICAL INNOVATION



Centrifugal vs Rotary Lobe

Centrifugal Pump

- Velocity Machine.
- Generates flow and head by the centrifugal action of a spinning impeller.



Rotary Lobe Pump


- Fixed Displacement Machine.
- Delivers a fixed flow rate per revolution of the rotors.



ENGINEERING PHARMACEUTICAL INNOVATION



Centrifugal Pump Operation

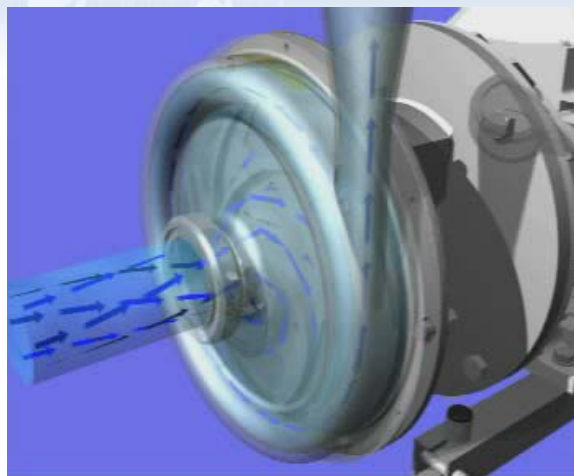

ENGINEERING PHARMACEUTICAL INNOVATION



Rotary Lobe Pump Operation



ENGINEERING PHARMACEUTICAL INNOVATION



Centrifugal vs Rotary Lobe

Centrifugal Pump

- Low Viscosity Fluid.
- Generally up to 500 Cps.
- Can go to 1500 Cps in special circumstances, but with large sacrifice in efficiency.
- Medium to high flows. Medium to low pressures.
- Both direct drive (1800 and 3600 rpm) and variable speed operation are common.
- Impeller diameter is usually trimmed to meet desired flow and head (pressure).

Rotary Lobe Pump

- High Viscosity Fluids.
- Generally 500 Cps and higher.
- Can use in low viscosity applications with sacrifice in efficiency.
- Low to medium flows. Medium to high pressures.
- Precise flow volume control, especially low flow rates.
- Generally operates at low speeds (up to 500 rpm) with a variable speed gear-motor.

ENGINEERING PHARMACEUTICAL INNOVATION



Control of Pumps

ENGINEERING PHARMACEUTICAL INNOVATION




Centrifugal Pumps:

A centrifugal pump will operate where the Pump Flow-Head Curve intersects with the System Resistance Curve.

Example:
DI Water
120 gpm
100 ft tdh (43 psi)
1.0 spgr

System Head:
25 ft Static
100 ft @ 120 gpm

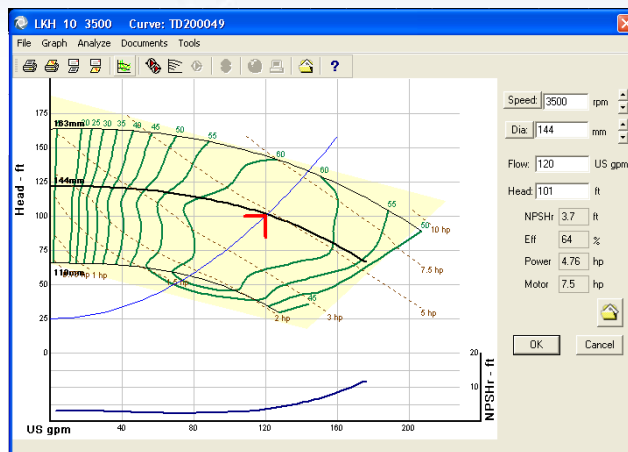

ENGINEERING PHARMACEUTICAL INNOVATION




Centrifugal Pumps:

Flow control – throttling the discharge with a control valve & flow meter.

Example:
DI Water
120 gpm
100 ft tdh (43 psi)
1.0 spgr

Pump:
3500 rpm
144mm impeller

Throttle:
80 gpm
55 gpm

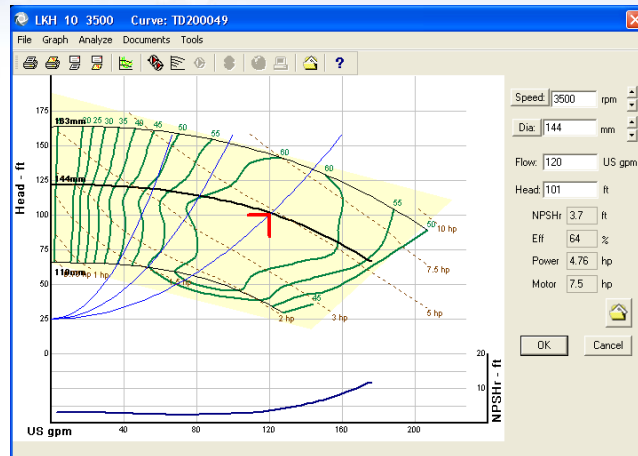


Centrifugal Pumps:

Flow control – variable speed drive (control motor rpm).

Example:
DI Water
120 gpm
100 ft tdh (43 psi)
1.0 spgr

Pump:
3500 rpm
144mm impeller

Speed Change:
80 gpm – 2550 rpm
55 gpm – 2000 rpm

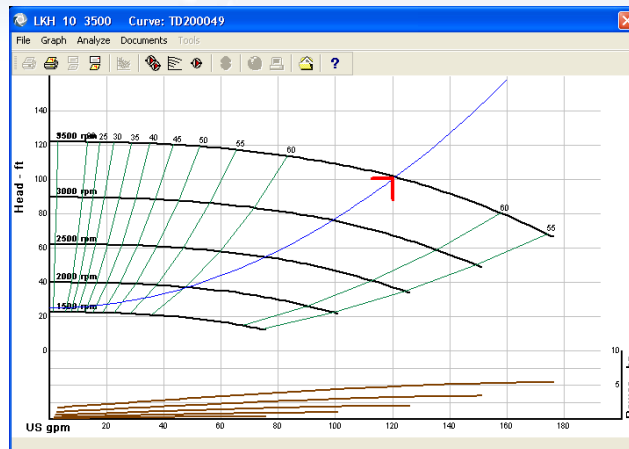

Rotary Lobe Pumps:

Flow control – variable speed drive (control motor rpm).

Example:

Fluid XYZ
30 gpm
100 psi
600 Cps
1.0 spgr

Rotary Lobe Pump

Model ABC
9.25 gal / 100 Rev.
0.0925 gal / Rev.

30 gpm = 324 rpm
25 gpm = 270 rpm
20 gpm = 216 rpm

Note: an RPM adjustment would be required for low viscosity fluids due to pump slip.

ENGINEERING PHARMACEUTICAL INNOVATION



Materials of Construction and Certifications

ENGINEERING PHARMACEUTICAL INNOVATION



Materials of Construction:

Application Dependent – For Example:

Pump Construction & Metallurgy

- Wetted Parts: 316L ss.
- Elastomers: EPDM, FPM, FDA Compliant, USP Class VI
- Finish: 32Ra, 20Ra, 15Ra, 10Ra, Electropolish, Passivation, etc.
- Mechanical Seal: Single or Double; Balanced vs Unbalanced; Face Materials; Flush or No Flush.
- Exterior: Painted or No Paint.
- Drainability, Cleanability, CIP-ability (more later).
- Standards: 3A, BPE, FDA, EHEDG, cGMP

Certifications & Tests

- Certifications: Material Traceability
- Performance Guarantees / Factory Tests: Q, H, E, P, N, etc.

ENGINEERING PHARMACEUTICAL INNOVATION



Mechanical Seals


ENGINEERING PHARMACEUTICAL INNOVATION



Mechanical Seals:

Hydraulically Balanced Seal Design:

- Important for pressure spikes.
- Important for pressure reversals.
- External seal (outside of process).

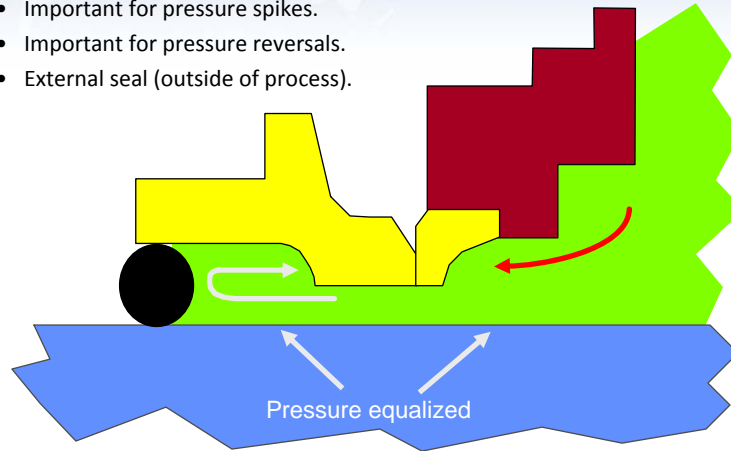

ENGINEERING PHARMACEUTICAL INNOVATION



Suction Condition

ENGINEERING PHARMACEUTICAL INNOVATION



Suction Conditions:

Huge Area for Potential Operational Problems:

Rule of Thumb:

- Fluid need to be able to get to pump suction, under sufficient pressure (no vaporization), as fast as the pump want to take it, or the pump will cavitate.

General Guidelines:

- Keep suction lines as short as practical. Short is good.
- Keep suction line at least one size larger than pump discharge. Bigger is better.
- Keep suction line velocity low. Especially with viscous fluids.
- Minimize pressure drop in suction line.
- Avoid air traps in suction line. Can cause loss of prime.

ENGINEERING PHARMACEUTICAL INNOVATION



Suction Conditions:

If You are Involved in Sizing & Selecting a Pump for a New Application:

Take Careful Note of the NPSH Required by the Pump:

- $NPSH = \text{Net Positive Suction Head}$.
- $NPSHR$ = Suction Head (pressure) required at the pump to prevent cavitation.
- Do an NPSH Available calculation for your specific piping arrangement.

NPSHA must be greater than NPSHR, or you will be Not Pumping So Hot !

What is Cavitation:

- Vaporization of fluid at eye of impeller.
- As the fluid & bubbles travels along the vane, the pressure increases.
- The collapse of those gas bubbles back to fluid is a high energy impact.
- Enough to remove metal from impeller. Deflect shafts. Rotor strike.

ENGINEERING PHARMACEUTICAL INNOVATION



Suction Conditions:

Cavitation Damage – Centrifugal Pump Impeller:

Signs of Cavitation:

- Centrifugal Pump – Sounds like pumping marbles, vibration, noise, damage.
- Rotary Lobe Pump – Banging noise (rotors hitting case), vibration.



ENGINEERING PHARMACEUTICAL INNOVATION



Pump Efficiency and Power Consumption

ENGINEERING PHARMACEUTICAL INNOVATION



Efficiency & Power consumption:

Positive Displacement Pumps

- Generally highly efficient. Typical 90 – 99%.
- Exception is with low viscosity fluids.


Centrifugal Pumps

- Generally range of efficiency – from < 10% (if running at very low flows) to 60 – 70%.
- Lots of room for careful review.
- The newer designs are generally 25% to 50% more efficient than the older pump technology.
- Reasons: Impeller geometry, improved manufacturing methods, tighter clearances, less slip (recirculation).

ENGINEERING PHARMACEUTICAL INNOVATION




Typical Centrifugal Pump:

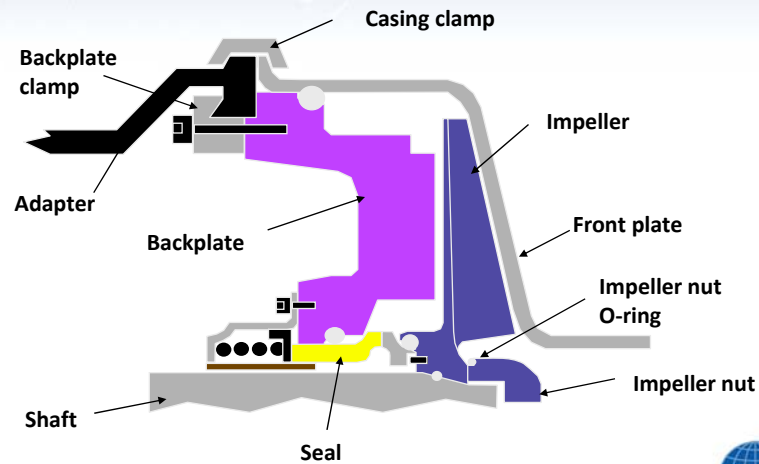

ENGINEERING PHARMACEUTICAL INNOVATION




Centrifugal Pump – Impeller Clearance:

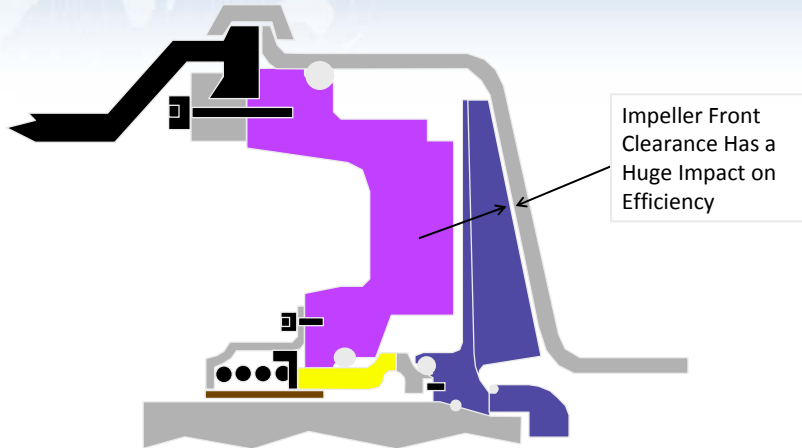

ENGINEERING PHARMACEUTICAL INNOVATION



Large Clearance = Recirculation (slip):

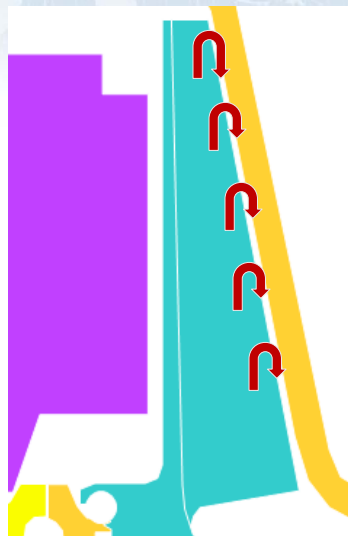

The greater the clearance between the impeller vanes and the pump volute, the higher the recirculation "slip" as more fluid simply churns within the pump, and the lower the resulting pump efficiency.

ENGINEERING PHARMACEUTICAL INNOVATION



Energy Savings:

Centrifugal Pumps – is it worth a look ?

Example: Centrifugal Pump - 100 gpm @ 100 ft tdh

Typical “old” Series	55% Eff	5.18 hp
High Eff. Pump	67% Eff	3.73 hp 28 % power reduction

At \$0.10/Kwh and 24/7 operation Power cost is \$4380/yr (High Eff)
Power cost is \$6084/yr (Old Series)
Savings = \$1704/yr.

Consider upgrading existing pumps too ?

Rule of thumb: 1 HP power reduction @ \$0.10/Kwh = \$1100/Yr. (cont. op).

ENGINEERING PHARMACEUTICAL INNOVATION




Cleanability
Drainability
CIP-ability

ENGINEERING PHARMACEUTICAL INNOVATION



Centrifugal Pumps

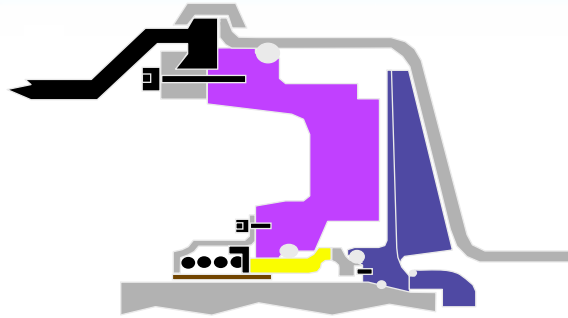


ENGINEERING PHARMACEUTICAL INNOVATION



Rotary Lobe Pumps

Drainability & Cleanability

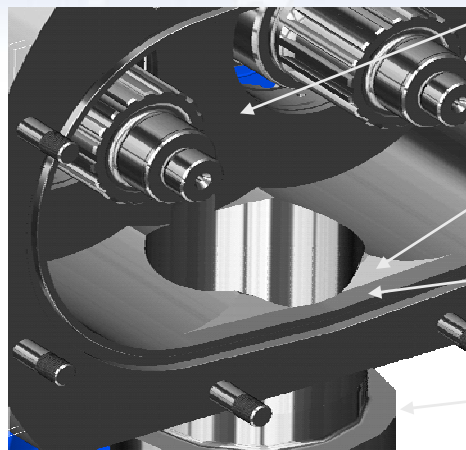

Smooth, Flat, Flush front and back surfaces.

Cusps set back

3° Deg fall to outlet

Vertical Porting


ENGINEERING PHARMACEUTICAL INNOVATION



Rotary Lobe Pumps

Drainability & Cleanability

Zero Crevice design
Front Cover Gasket

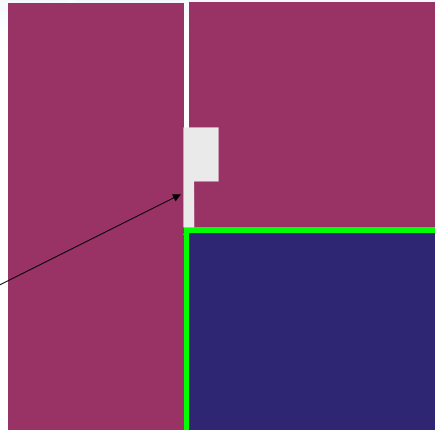

ENGINEERING PHARMACEUTICAL INNOVATION



Cleanability, Drainability, CIP-ability:

Centrifugal Pumps

- 3A, FDA, EHEDG certifications.
- Internal finish.
- Low point drains.
- Design & position of o-rings & gaskets.

Positive Displacement Pumps

- Look beyond the 3A certification.
- Vertical porting.
- Drainable cusps.
- Smooth, flat, flush rotor case back & front.
- Mechanical seal position for cleanability.

ENGINEERING PHARMACEUTICAL INNOVATION





Centrifugal and Positive Displacement Pumps For Hygienic Applications

Thanks for your attention.

Questions.

Presented by:

Allen E. LeBoeuf, Product Manager

Oliver M. Dean, Inc. – Worcester, MA

ENGINEERING PHARMACEUTICAL INNOVATION

